




WARNING! DEATH or SERIOUS INJURY can occur if manufacturer's instructions are not followed correctly

Product Instructions

E-Z-ON Portable Ankle Restraints : Item #100AR

- * **Use** for added security when extra restraint is needed.
- * **Three McDonald buckle** closures for quick securement and easy one-pull adjustment
- * **Deter** unwanted leg movements that can endanger others.
- * **Considered** a less restrictive restraint.
- * **Wrap** around each ankle and secure push-buckle closures.
- * **Position** corresponding McDonald buckle on inside of ankle and secure to each other.
- * **Adjust** by pulling on excess webbing to secure restraints snug around ankle.
- * **Back thread** webbing through outside push-buckle slot to lock into place.
- * **Follow** these instructions for proper use and care.

Note: Sold as shown

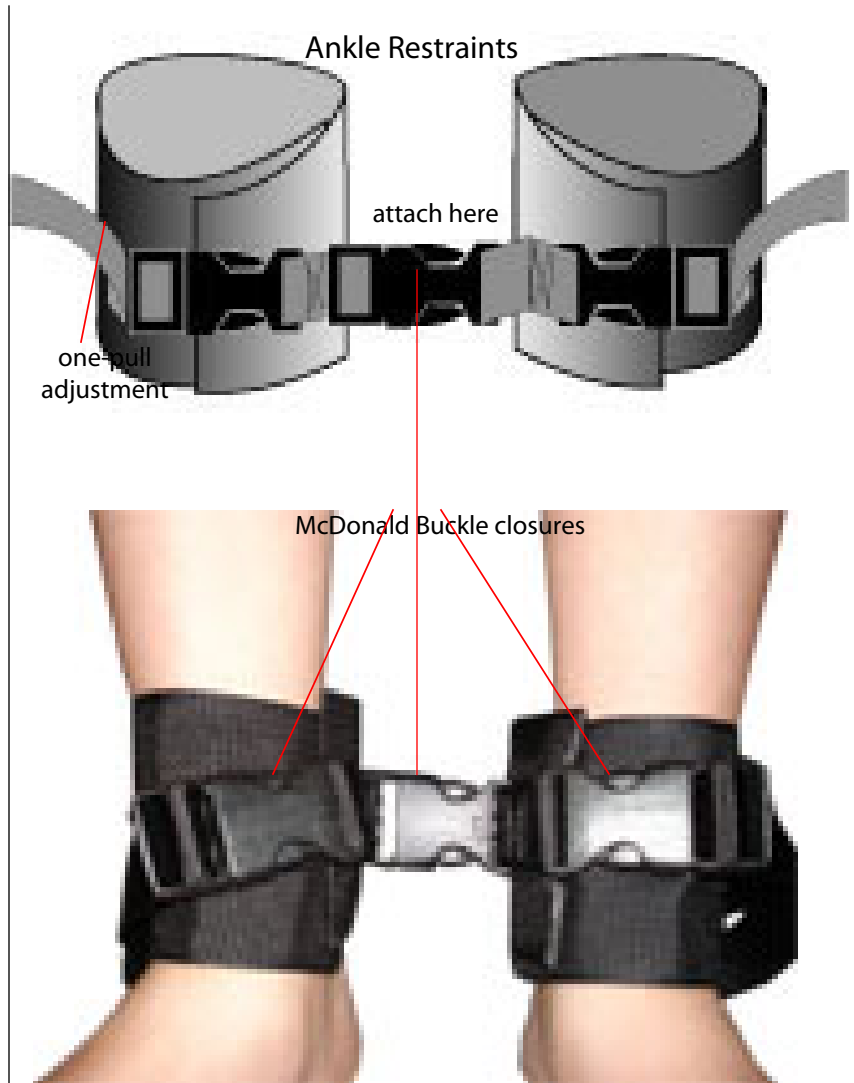

E-Z-ON PRODUCTS INC. OF FLORIDA
605 Commerce Way West
Jupiter, FL 33458
800-323-6598 FAX 561-747-8779